
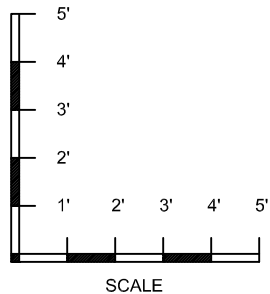


MOBILE X-RAY BARRIER.
CUSTOMER SUPPLIED,
CUSTOMER INSTALLED.

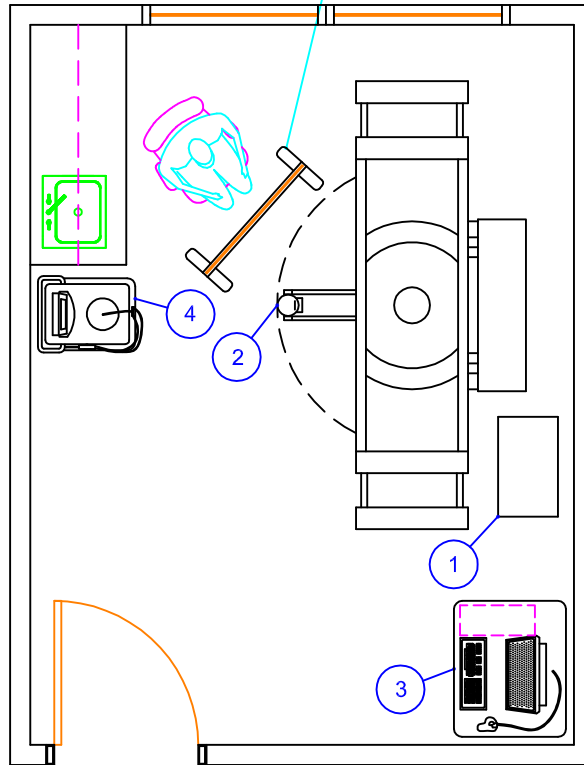

TYPICAL ROOM AND EQUIPMENT LAYOUT

GENERAL NOTES

• ALL SPECIFICATIONS AND SPACE REQUIREMENTS ARE SHOWN GENERICALLY. CONTACT GAMMA TECH, INC. FOR SPECIFIC REQUIREMENTS PERTAINING TO OPTIONAL EQUIPMENT COMPONENTS AND ADDITIONAL ROOM PREPARATION SPECIFICATIONS.

VOLTAGE: -10%, +10% 200,208,220,230,240VAC
 CURRENT: (INPUT) 4A/IDLE,35A/5 SECONDS MAX @ 220VAC
 FREQUENCY: ± .5 Hz 50/60 Hz
 MAX POWER: INPUT 5KW
 MAX KVA: (OUTPUT) 2.38KVA (70mA@34kV)
 LINE TYPE: DEDICATED HARDWIRED FOR 30AMP SERVICE
 CIRCUIT BREAKER: DUAL LINE, DEDICATED, 15 AMP ON GENERATOR

EQUIPMENT SCHEDULE

EQUIPMENT IDENTIFICATION NUMBER

		BTU / HOUR		WEIGHT (LBS.)		
		DESCRIPTION		DIMENSIONS (IN.)		
		W	D	x	H	
1	1100	LORAD MULTICARE PLATINUM GENERATOR	25.00	15.00	20.00	190
2		LORAD MULTICARE PLATINUM BREAST BIOPSY TABLE	82.50	66.00	57.00	600
3		LORAD MULTICARE PLATINUM DSM VIEWING AND PROCESSING CONSOLE	34.00	28.00	45.00	257
4		ETHICON MAMMOTOME SYSTEM CONTROL MODULE W/ SMARTVAC™	18.00	24.00	48.00	

GAMMA TECH, INC



129 Wild Berry Lane
 Hampstead, NC 28443
 Ph (910) 270-9787
 Fax (910) 270-9767
 www.gammatech.net

© PLANNING DIVISION
 FILE NAME:
 LORAD MULTICARE PM

REVISION DATE:
 07.17.2008

LORAD

MULTICARE PLATINUM
 BIOPSY SUITE