



# EQUIPMENT SCHEDULE

EQUIPMENT IDENTIFICATION NUMBER						
		BTU / HOUR	WEIGHT (LBS.)			
		DESCRIPTION	DIMENSIONS (IN.) W x D x H			
1	12	SHIMADZU UD150B-30 X-RAY CONTROL CONSOLE MOUNTED ON COUNTER	13.20	11.40	3.15	12
2	3095	SHIMADZU UD150B-30 X-RAY CONTROL CABINET	31.50	15.80	68.90	551
3	1800	AMERICAN TABLE INC. ATM 6000 CV (TABLE LENGTH INCLUDES EXTENSIONS)	156.00	27.00	VAR.	990
4	3817	SHIMADZU MH-200 WITH 2 METER EXTENDED RAILS CEILING SUSPENDED C-ARM SYSTEM	295.25 (MODIFIED 236.25)	114.32	112.25	2640
5	794	SHIMADZU MH-200 CONTROL CABINET	23.62	17.72	71.90	465
6	342	H.M.S. ELECTRONICS HS3-MPX HIGH SPEED STARTER	20.00	9.00	17.00	70
7		IMAGE DIAGNOSTIC, INC. MODEL ID1620RF-34.5 CEILING SUSPENDED DUAL MONITOR SYSTEM	160.00 (MODIFIED 134.00)	17.12	55.00	200
8	680	SHIMADZU 21" IMAGE DISPLAY MONITOR (X2) (CEILING SUSPENDED)	19.30	21.60	19.30	55
9	410	SHIMADZU 17" IMAGE DISPLAY MONITOR MOUNTED ON THE COUNTER	16.50	16.50	16.10	44
10	5120	SHIMADZU DIGITEX PRO MULTI DIGITAL COMPUTER MOUNTED UNDER THE COUNTER	8.70	22.00	23.60	77
11	N/A	SHIMADZU DIGITEX PRO MULTI POWER DISTRIBUTION BLOCK MOUNTED UNDER THE COUNTER	15.90	3.40	6.50	5
12	800	SHIMADZU DIGITEX PRO MULTI CAMERA CONTROL MOUNTED ON SHELF UNDER THE COUNTER	14.20	17.80	5.90	25
13	6800	SHIMADZU SERIES ST, MODEL 1020 UPS (2 KVA) MOUNTED UNDER THE COUNTER	7.00	18.00	13.00	81
14	1300	SHIMADZU PDCU-SRF 130KVA POWER CONDITIONER	21.70	20.50	44.50	730
15	1200	CODONICS Z-200 NP-1660M MEDICAL PRINTER	17.00	21.00	12.00	69
16	56	SHIMADZU CD-R IMAGE ARCHIVE SYSTEM MOUNTED ON COUNTER	6.70	13.00	2.70	7

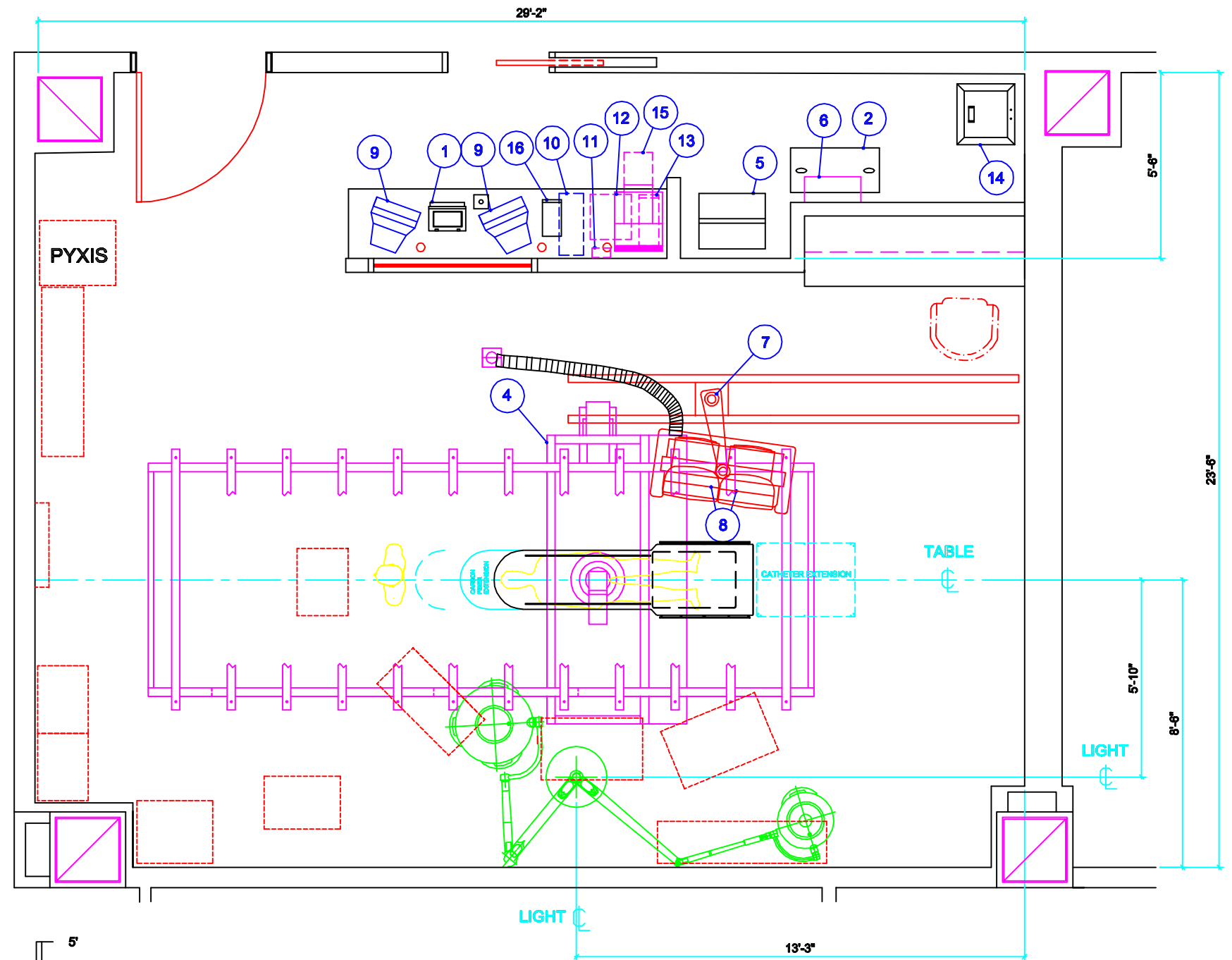

**PROPOSED EQUIPMENT LAYOUT**  
(DASHED LINES DEPICT TABLE TOP TRAVEL)  
(EQUIPMENT SHOWN IN IMAGING PROCEDURE POSITION)

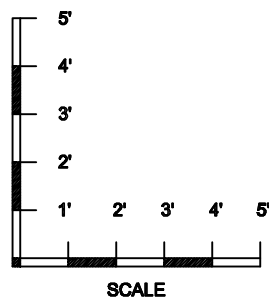

RECOMMENDED POWER TO THE MAIN DISCONNECT: 480VAC, 60Hz, THREE PHASE, 130KVA (DELTA OR WYE)  
 MAIN CIRCUIT BREAKER: 110AMP, 3-POLE, THERMOMAGNETIC  
 MAXIMUM MOMENTARY: 180 AMPS  
 X-RAY GENERATOR: SHIMADZU MH-200 SYSTEM WITH UD150B-30 AND PDCU-SRF

**SHIMADZU**

**MH-200 DIGITAL  
ENDOVASCULAR ROOM**

( TOTAL SQUARE FOOTAGE = 686 )  
( RECOMMENDED CEILING HEIGHT = 10'-0" )

PLANNING BY:



**GAMMA TECH, INC**

129 Wild Berry Lane  
Hampstead, NC 28443  
Ph (910) 270-9787 Fax (910) 270-9767  
www.gammatech.net